
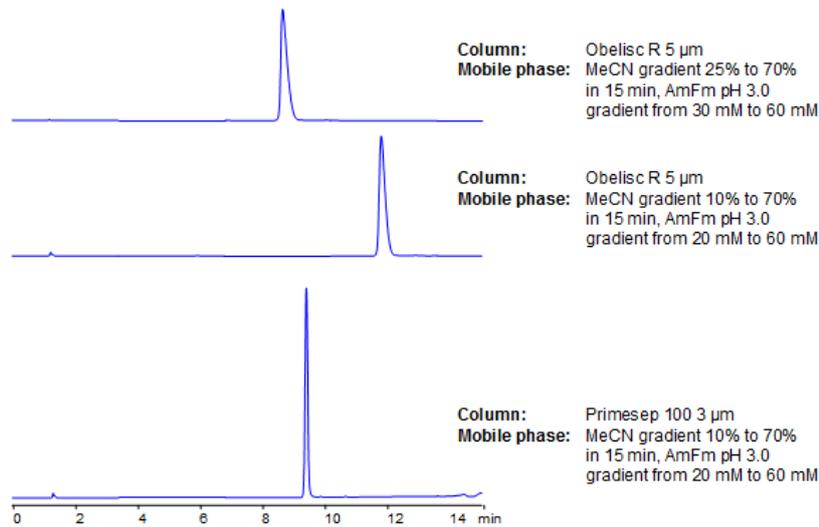


## 1-Naphthaleneacetic acid (NAA) Separation on Primesep and Obelisc Mixed-Mode Columns

Size: 2.1 x 150 mm  
 Flow: 0.4 mL/min  
 Detection: UV 270 nm

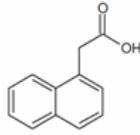

SIELC has developed the Obelisc™ columns, which are mixed-mode and utilize Liquid Separation Cell technology (LiSC™). These cost-effective columns are the first of their kind to be commercially available and can replace multiple HPLC columns, including reversed-phase (RP), AQ-type reversed-phase, polar-embedded group RP columns, normal-phase, cation-exchange, anion-exchange, ion-exclusion, and HILIC (Hydrophilic Interaction Liquid Chromatography) columns. By controlling just three orthogonal method parameters - buffer concentration, buffer pH, and organic modifier concentration - users can adjust the column properties with pinpoint precision to separate complex mixtures.

### Method Parameters

<b>Column</b>	Primesep 100, 2.1x150 mm, 3 µm, 100 Å
<b>Mobile Phase</b>	Gradient MeCN – 10-70%, 15 min
<b>Buffer</b>	Gradient AmAc pH 3.0- 20-60 mM, 15 min
<b>Flow Rate</b>	0.4 mL/min
<b>Detection</b>	UV, 270 nm

Quelle: <https://sielc.com/compound-1-naphthalene>