


Analysis of 2,4,5-Trichlorophenoxyacetic acid on Primesep B2 Column

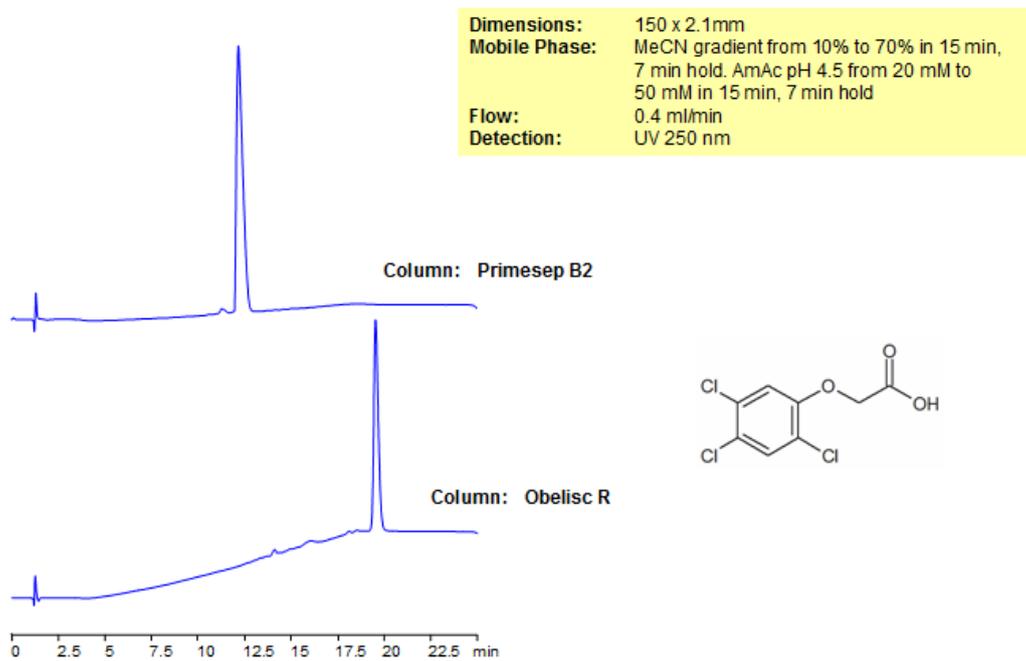

Method Parameters

Column	Primesep B2, 2.1x150 mm, 5 µm, 100 Å
Mobile Phase	Gradient MeCN – 10-70%, 15 min, 7 min hold
Buffer	Gradient AmAc pH 4.5- 20-50 mM, 15 min, 7 min hold
Flow Rate	0.4 mL/min
Detection	UV, 250 nm

Quelle:

<https://sielc.com/Application-Analysis-of-2-4-5-Trichlorophenoxyacetic-acid-on-Primesep-and-Obelisc-HPLC-Columns>