


Benfluralin Analysis on Obelisc R HPLC Column

Column:	Obelisc R
Dimensions:	150 x 2.1mm
Mobile Phase:	MeCN gradient from 10% to 70% in 15 min, 7 min hold. AmAc pH 4.5 from 20 mM to 50 mM in 15 min, 7 min hold
Flow:	0.4 ml/min
Detection:	UV 250 nm

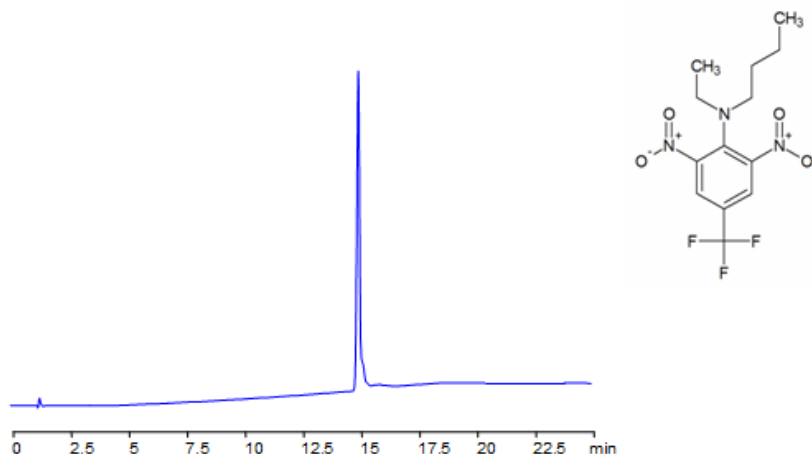

SIELC has developed the Obelisc™ columns, which are mixed-mode and utilize Liquid Separation Cell technology (LiSC™). These cost-effective columns are the first of their kind to be commercially available and can replace multiple HPLC columns, including reversed-phase (RP), AQ-type reversed-phase, polar-embedded group RP columns, normal-phase, cation-exchange, anion-exchange, ion-exclusion, and HILIC (Hydrophilic Interaction Liquid Chromatography) columns. By controlling just three orthogonal method parameters - buffer concentration, buffer pH, and organic modifier concentration - users can adjust the column properties with pinpoint precision to separate complex mixtures.

Method Parameters

Column	Obelisc R, 2.1×150 mm, 5 µm, 100 Å
Mobile Phase	Gradient MeCN – 10-70%, 15 min, 7 min hold
Buffer	Gradient AmAc pH 4.5- 20-50 mM, 15 min, 7 min hold
Flow Rate	0.4 mL/min
Detection	UV, 250 nm

Quelle: <https://sielc.com/Application-Benfluralin-Analysis-on-Obelisc-R-HPLC-Column>