


HPLC Method for Analysis of 4-Amino-2-nitrophenol on Primesep 100 Column

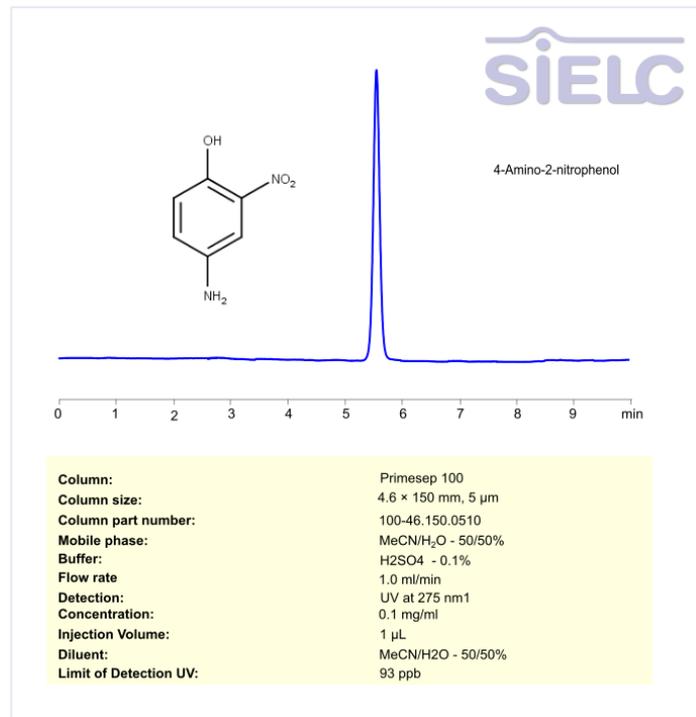

High Performance Liquid Chromatography (HPLC) Method for Analysis of 4-Amino-2-nitrophenol

4-Amino-2-nitrophenol is an organic compound with the molecular formula C₆H₆N₂O₃.

Properties: Appearance: Typically a yellow to orange crystalline solid

Molecular weight: ~154.12 g/mol

Solubility: Soluble in water and polar organic solvents

Uses: Intermediate in the manufacture of dyes, pigments, and pharmaceuticals

4-Amino-2-nitrophenol can be retained and analyzed using the Primesep 100 stationary phase column. The analysis utilizes an isocratic method with a simple mobile phase consisting of water, acetonitrile (MeCN), and sulfuric acid. Detection is performed using UV at 275 nm.

*LOD was determined for this combination of instrument, method, and analyte, and it can vary from one laboratory to another even when the same general type of analysis is being performed.

Method Parameters

Column	Primesep 100, 4.6 x 150 mm, 5 µm, 100 Å, dual ended
Mobile Phase	MeCN – 50%
Buffer	H2SO4 – 0.1%
Flow Rate	1.0 mL/min
Detection	UV 275 nm
Limit of Detection	93 ppb

Quelle: <https://sielc.com/hplc-ms-method-for-analysis-4-amino-2-nitrophenol>