


HPLC Method for Analysis of Iprodione on Newcrom B Column on Alltesta™

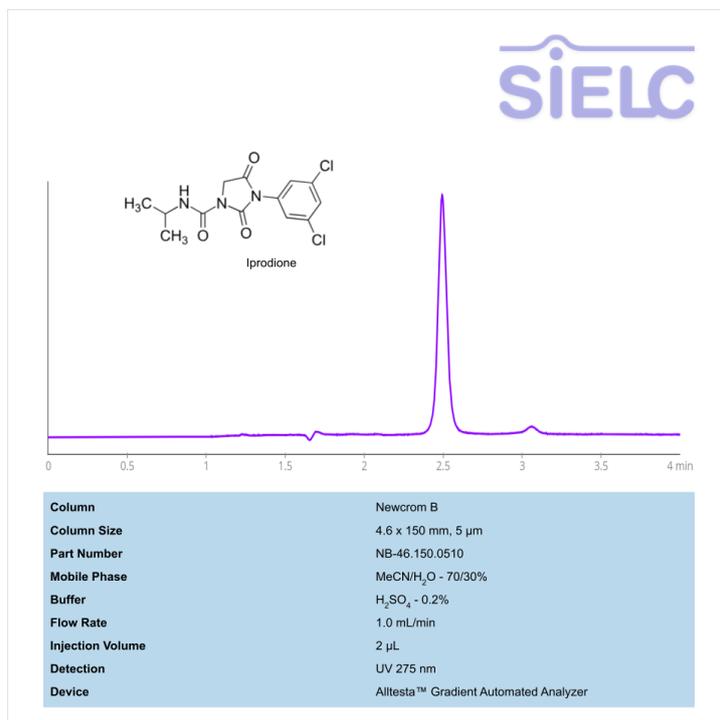

High Performance Liquid Chromatography (HPLC) Method for Analysis of Iprodione

Iprodione is a hydantoin fungicide and nematicide with the chemical formula C₁₃H₁₃Cl₂N₃O₃. It is used on a wide variety of crops that are affected by fungal diseases such as botrytis bunch rot, brown rot, and Sclerotinia. Approval of using it is different across the world.

Iprodione can be retained and analyzed using the Newcrom B stationary phase column. The analysis utilizes an isocratic method with a simple mobile phase consisting of water and acetonitrile (MeCN) with phosphoric acid as a buffer. Detection is performed using UV.

Method Parameters

Column	Newcrom B, 4.6 x 150 mm, 5 µm, 100 Å, dual ended
Mobile Phase	MeCN – 70%
Buffer	Sulfuric Acid
Flow Rate	1.0 mL/min
Detection	UV 275 nm

Quelle: <https://sielc.com/hplc-method-for-analysis-of-iprodione>