


## HPLC Method for Analysis of Lipitor

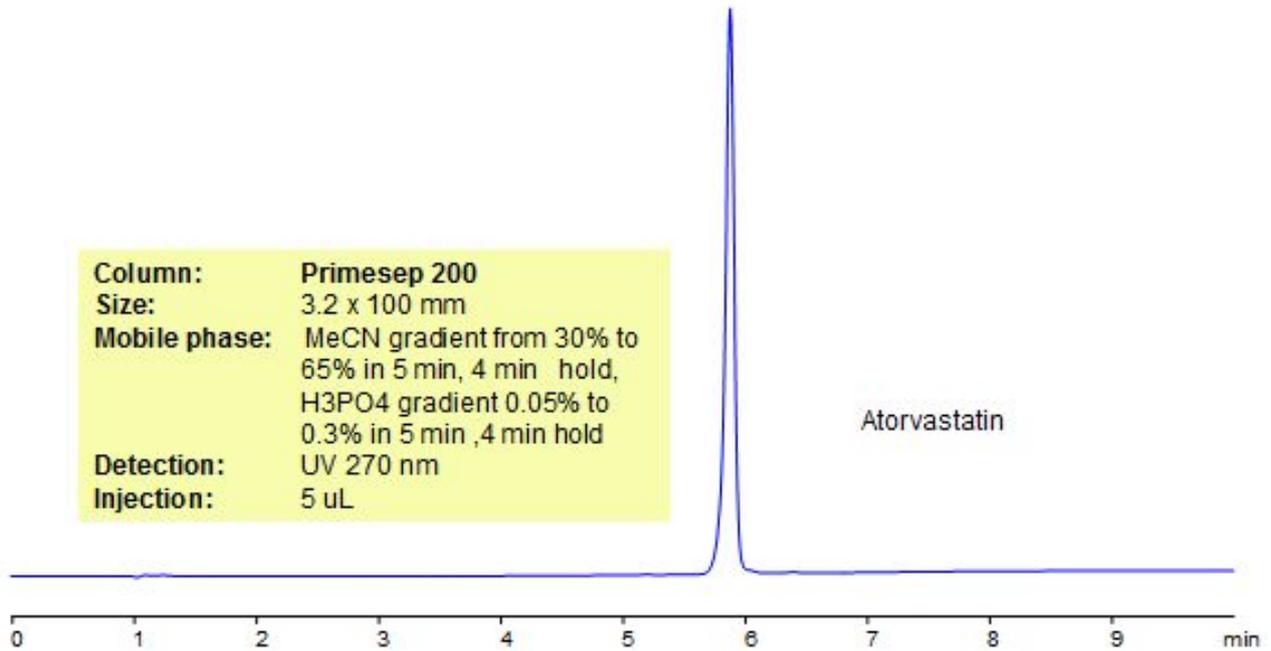

### Method Parameters

<b>Column</b>	Primesep 200, 3.2x100 mm, 5 µm, 100 Å
<b>Mobile Phase</b>	Gradient MeCN – 30-65%, 5 min, 4 min hold
<b>Buffer</b>	Gradient H3PO4 0.05 – 0.3% , 5 min, 4 min hold
<b>Flow Rate</b>	1.0 mL/min
<b>Detection</b>	UV, 270 nm

Quelle: <https://sielc.com/HPLC-Method-Analysis-Lipitor>