


HPLC Method for Analysis of Nirmatrelvir in PAXLOVID™ (Nirmatrelvir Tablets) on Primesep P Column

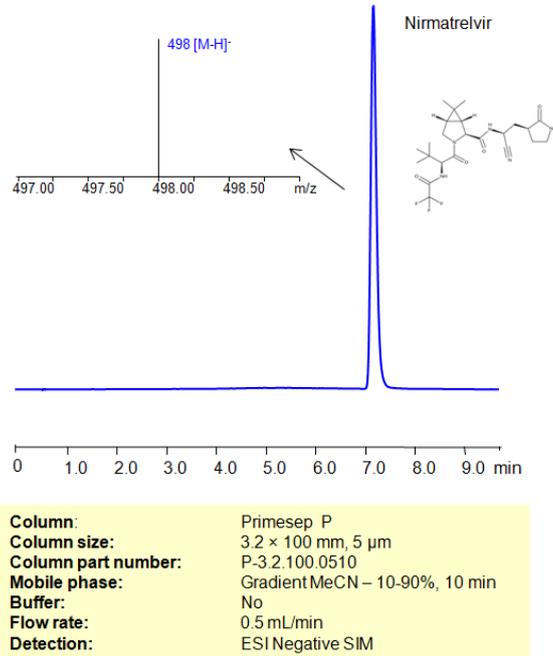

Separation type: Liquid Chromatography Mixed-mode

High Performance Liquid Chromatography (HPLC) Method for Analysis of Nirmatrelvir

Method Parameters

Column	Promesep P, 3.2×100 mm, 100 Å
Mobile Phase	Gradient MeCN – 10-90%
Buffer	No
Flow Rate	0.5 mL/min
Detection	UV 210 nm

Quelle: <https://sielc.com/hplc-method-for-analysis-of-nirmatrelvir-in-paxlovid-nirmatrelvir-tablets>