
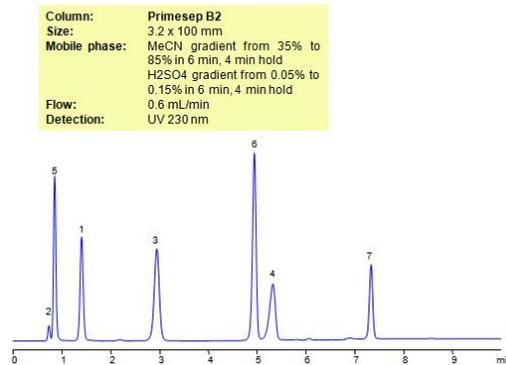


HPLC Method for Analysis of Pesticides: Thiamphenicol, Aniline, Simazine, Alachlor, Thiabendazole, Diazinon, Trifluralin on Primesep B2 Column

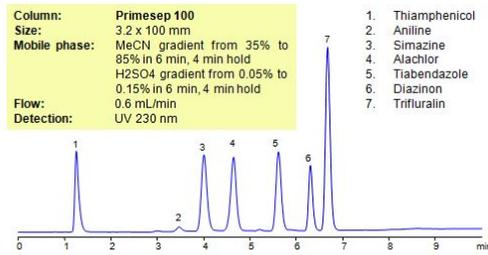

Method Parameters

Column	Primesep B2, 3.2x100 mm, 5 µm, 100 Å
Mobile Phase	Gradient MeCN – 35-85%, 6 min , 4 min hold
Buffer	Gradient H ₂ SO ₄ – 0.05- 0.15%, 6 min, 4 min hold
Flow Rate	0.6 mL/min
Detection	UV, 230 nm

Quelle: <https://sielc.com/separation-of-pesticides>