


## HPLC Method for Analysis of Vicodin (Acetaminophen & Hydrocodone)

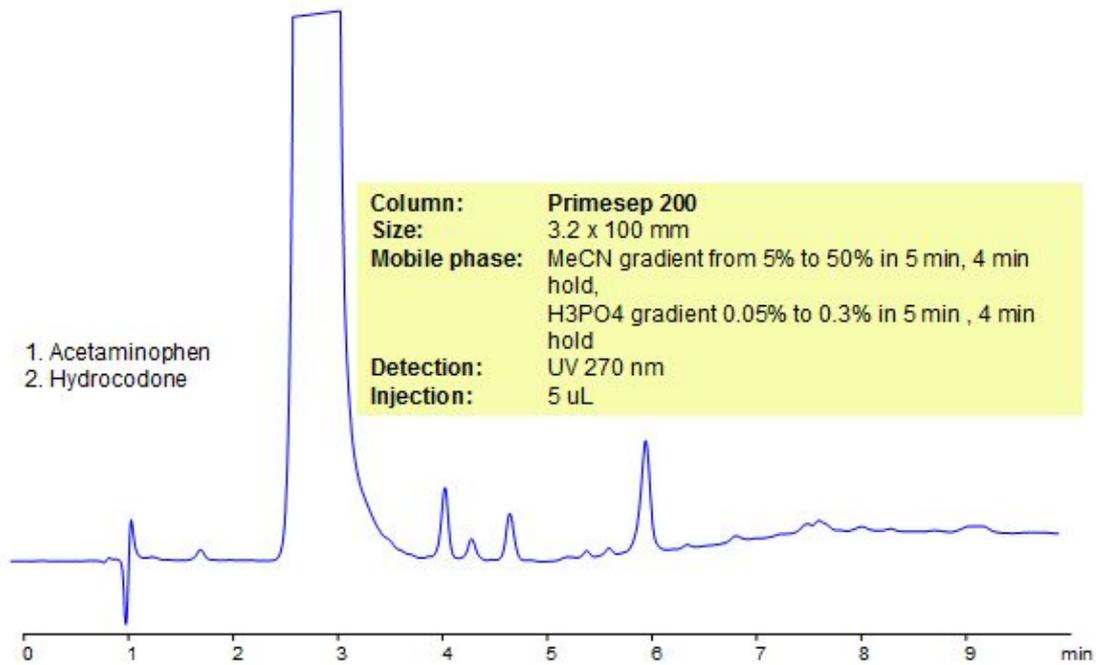

### Method Parameters

<b>Column</b>	Primesep 200, 3.2x100 mm, 5 $\mu$ m, 100 Å
<b>Mobile Phase</b>	Gradient MeCN – 5-50%
<b>Buffer</b>	Gradient H3PO4 – 0.05-0.3%, 5 min, 4 min hold
<b>Flow Rate</b>	0.5 mL/min
<b>Detection</b>	UV, 270 nm

Quelle: <https://sielc.com/HPLC-Method-Analysis-Vicodin>