


## HPLC Method for Separation of Hydrophobic, Cationic, Nonionic and Anionic Surfactants on Newcrom BH Column

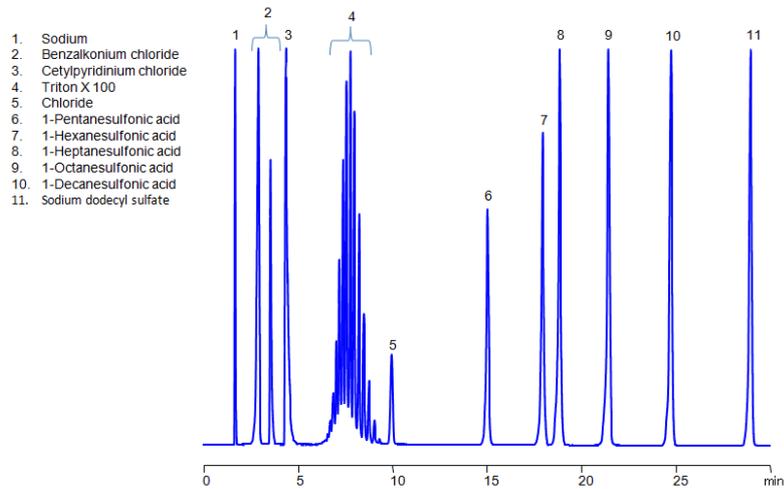

<b>Column:</b>	Newcrom BH
<b>Column size:</b>	4.6 × 250 mm, 5 µm, 100Å
<b>Column part number:</b>	NBH-46.250.0510
<b>Mobile phase:</b>	Gradient MeCN – 40-80%, 30 min
<b>Buffer:</b>	Ammonium formate pH 3.0 – 20 mM
<b>Flow rate:</b>	1.0 mL/min
<b>Detection:</b>	ELSD, 50C

Separation type: Liquid Chromatography Mixed-mode

Surfactants, also known as surface-active agents, are compounds that lower the surface tension (or interfacial tension) between two liquids or between a liquid and a solid. Surfactants may act as detergents, wetting agents, emulsifiers, foaming agents, or dispersants.

They are often classified according to the charge of the polar head group:

**Anionic Surfactants:** These surfactants have a negative charge on their polar head group. Common examples include soap, sodium laureth sulfate, and sodium lauryl sulfate. They are commonly used in detergents and shampoos due to their ability to emulsify oils and hold dirt in suspension, so it can be rinsed away.

**Cationic Surfactants:** These surfactants have a positive charge on their polar head group. Examples include cetyltrimethylammonium bromide (CTAB) and benzalkonium chloride. These are often used as antiseptics and can also be found in hair conditioners because they reduce static cling.

**Nonionic Surfactants:** These surfactants have no charge on their polar head group. Examples include alcohol ethoxylates, nonylphenol ethoxylates, and polysorbates. Nonionic surfactants are often used in laundry and dishwasher detergents.

All compounds can be retained, separated, and analyzed using a reverse-phase Newcrom BH, 4.6 x 250 mm, 5 µm, 100 Å, dual ended column. The mobile phase for this method consists of water, acetonitrile (MeCN), and Ammonium formate, which serves as a buffer. This analytical method can be detected with an Evaporative Light Scattering Detector (ELSD) or any other evaporative detection method (CAD,

ESI-MS).

### Method Parameters

<b>Column</b>	Newcrom BH, 4.6 x 250 mm, 5 µm, 100 Å, dual ended
<b>Mobile Phase</b>	Gradient MeCN -40-80%, 30 min
<b>Buffer</b>	Ammonium formate pH 3.0 – 20 mM
<b>Flow Rate</b>	1.0 mL/min
<b>Detection</b>	ELSD, 50C

Quelle: <https://sielc.com/hplc-separation-surfactants>