


HPLC Separation of Amino Acids in HILIC and Cation-Exchange Modes

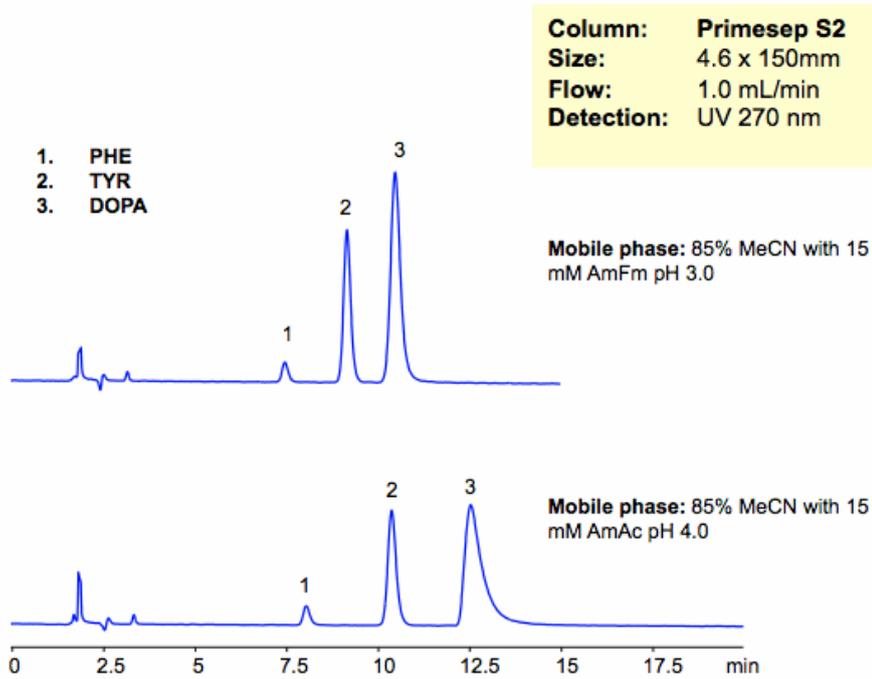

Amino acids can be retained and separated on a Primesep S2 column in HILIC/cation-exchange mode.

Method Parameters

Column	Primesep S2, 4.6x150 mm, 5 µm, 100 Å
Mobile Phase	MeCN/H ₂ O
Buffer	AmFm
Flow Rate	1.0 mL/min
Detection	UV, 270 nm

Quelle: <https://sielc.com/Application-HPLC-Separation-of-Amino-Acids-in-HILIC-and-Cation-Exchange-Modes>