


HPLC Separation of Gentamicins

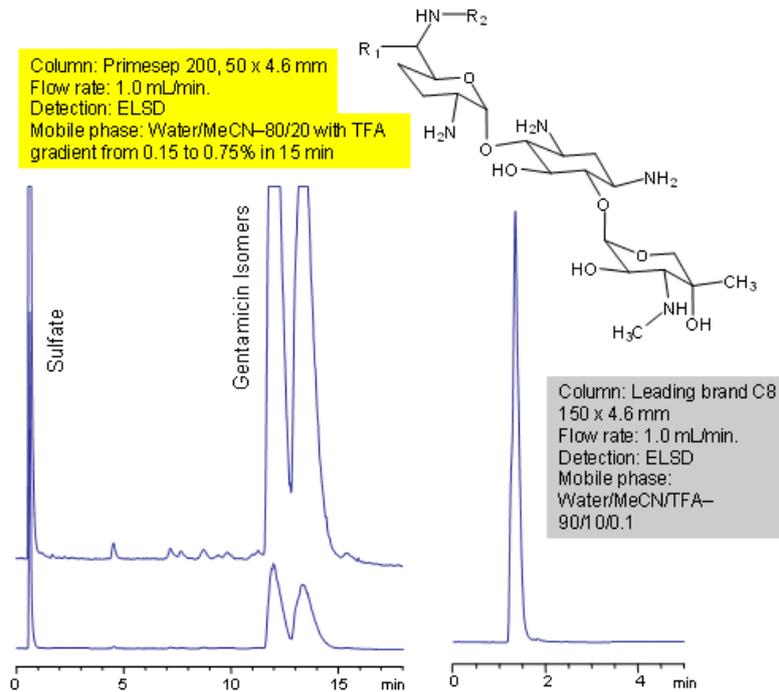

Primesep 200 separates the very hydrophilic gentamicin isomers. Gentamicin contains multiple primary and secondary amines as well as multiple hydroxyl groups. Ion-exchange and other polar interactions provide retention for a compound that elutes in the void of a traditional C8 column. The separation uses a simple mobile phase of water, acetonitrile (MeCN, ACN), and trifluoroacetic acid (TFA) with evaporative light scattering detection (ELSD).

Method Parameters

Detection	ELSD Detection
-----------	----------------

Quelle: <https://sielc.com/Application-HPLC-Separation-of-Gentamicins>