

HPLC Separation of Nitroprusside and its Metabolites on Primesep B2

1. Nitroprusside
2. Potassiumferrocyanide III
3. Potassiumferrocyanide II
4. Potassium thiocyanate

Column: Primesep B2

Size: 3.2 x 50 mm


Mobile phase:

Gradient from A to B in 5 min + 4 min hold B

Solvent A: 7 mM AmAc pH 4 with 25% MeCN

Solvent B: 7 mM AmAc pH 9 with 25% MeCN

Flow rate: 0.5 mL/min

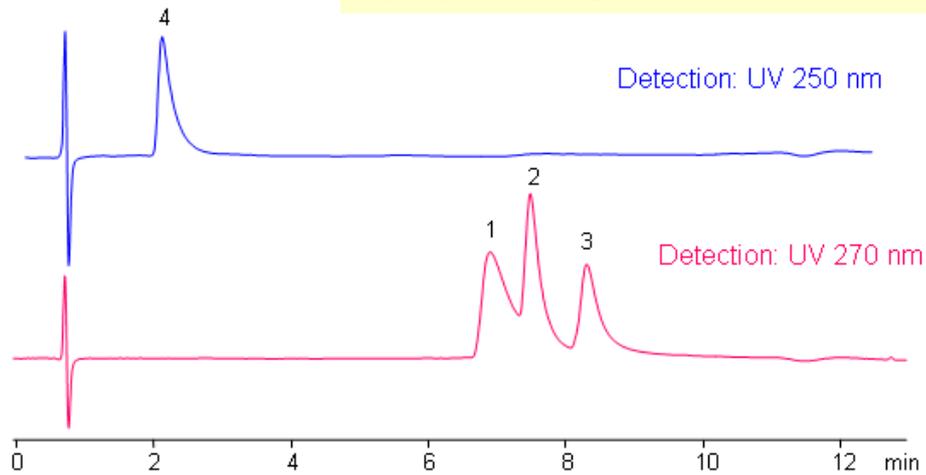

Method Parameters

Detection	UV Detection
-----------	--------------

Quelle: <https://sielc.com/Application-HPLC-Separation-of-Nitroprusside-and-Its-Metabolites>