


HPLC Separation of Nucleosides and Deoxynucleosides

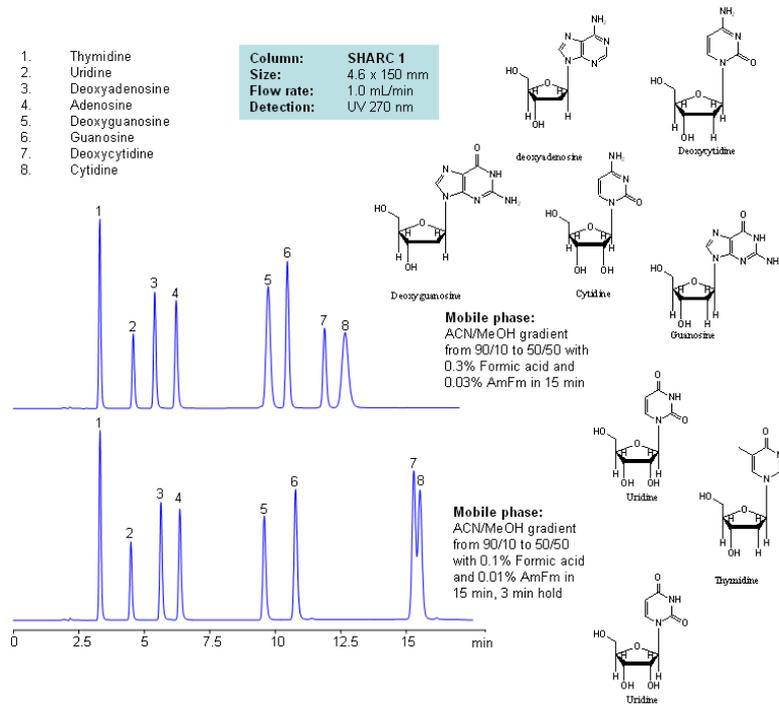

The SHARC™ family of innovative columns represents the first commercially available columns primarily utilizing separation based on hydrogen bonding. SHARC stands for Specific Hydrogen-bond Adsorption Resolution Column. Hydrogen bonding involves an interaction or attraction between a bound hydrogen atom and molecules containing electronegative atoms, such as oxygen, nitrogen, and fluorine.

Method Parameters

Column	Sharc 1, 4.6x150 mm, 5 µm, 100 Å
Mobile Phase	MeCN/MeOH
Buffer	AmFm, Formic acid
Flow Rate	1.0 mL/min
Detection	UV, 270 nm

Quelle: <https://sielc.com/Application-HPLC-Separation-of-Nucleosides-and-Deoxynucleosides>