


HPLC Separation of Succinylcholine and Choline Using a Primesep 200 Column

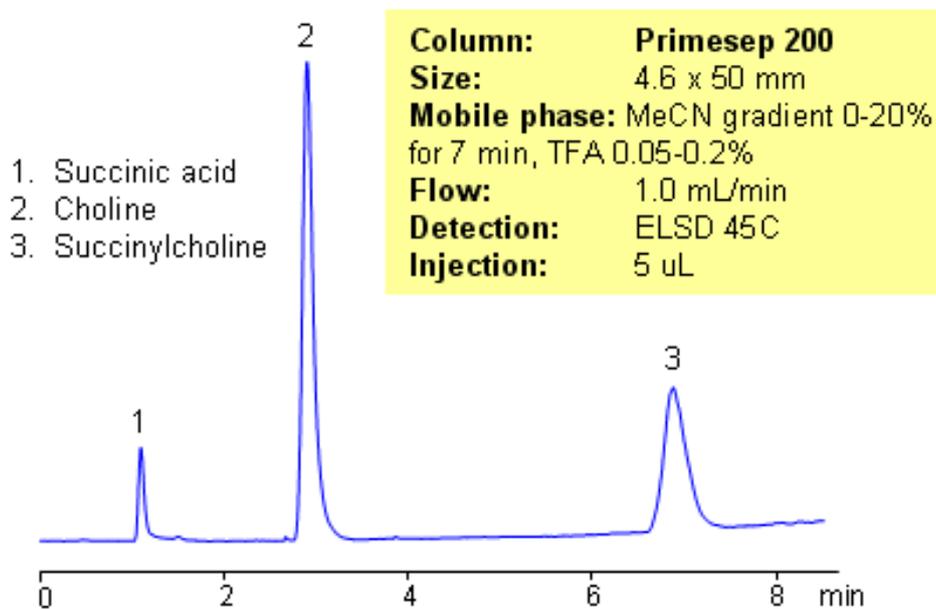

Method Parameters

Column	Primesep 200, 4.6x50 mm, 5 µm, 100 Å
Mobile Phase	Gradient MeCN0 – 20%, 7 min
Buffer	TFA – 0.05-0.20%
Flow Rate	1.0 mL/min
Detection	ELSD45C

Quelle: <https://sielc.com/Application-HPLC-Separation-of-Succinylcholine-and-Choline-Using-a-Primesep-200-Column>