


Reverse Elution Order in Separation of Alkali Cation

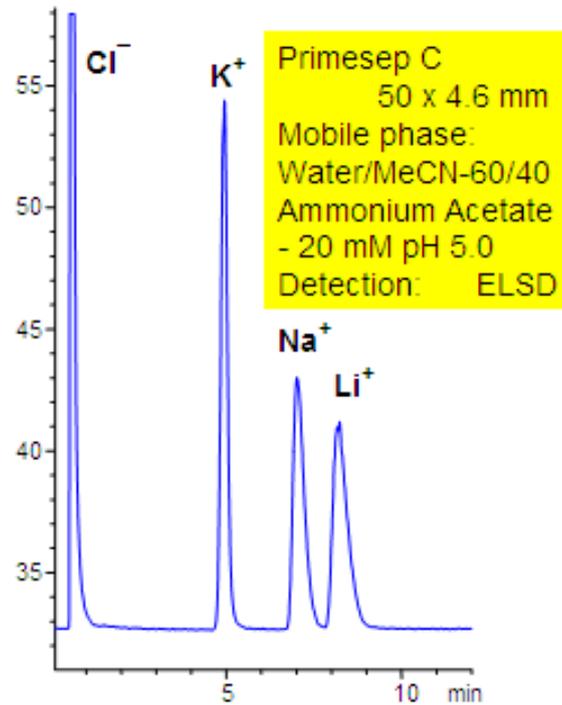

Primesep C separates a mixture of inorganic cations, potassium, sodium, and lithium, by a combination of cation exchange and complex formation. The peak order is reversed on the Primesep C column compared to the order normally obtained on columns that rely on cation exchange alone. The mobile phase mixture of water, acetonitrile (MeCN, ACN) and ammonium acetate is compatible with UV and evaporative light scattering detection (ELSD).

Method Parameters

Column	Primesep C, 4.6x50 mm, 5 µm, 100 Å
Mobile Phase	MeCN/H ₂ O
Buffer	AmAc
Flow Rate	1.0 mL/min
Detection	ELSD50C

Quelle: <https://sielc.com/Application-Reverse-Elution-Order-in-Separation-of-Alkali-Cation>