


HPLC Separation of Mono-, Di- and Trisaccharides on Amaze HD HILIC Mixed-Mode Column

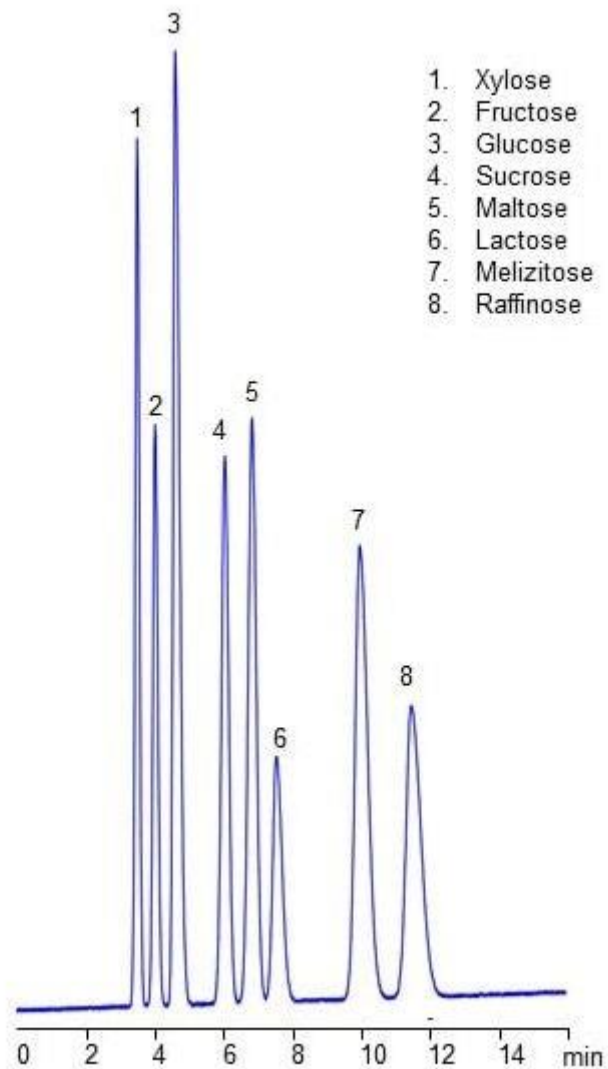

Column: Amaze HD

Separation Modes: HILIC, anion-exclusion

Column Dimensions: 4.6x150 mm, 3 μ m, 100A

Mobile Phase: ACN/water/AmAc pH 5

Detection: Corona CAD

Sample: 2 mg/ml

Injection: 2 μ L

Flow rate: 1 ml/min

Compounds: Xylose, Fructose, Glucose, Sucrose, Maltose, Lactose, Melizitose, Raffinose