


HPLC Method for Analysis of Vicodin (Acetaminophen & Hydrocodone)

<https://sielc.com/HPLC-Method-Analysis-Vicodin>

Chromatogram

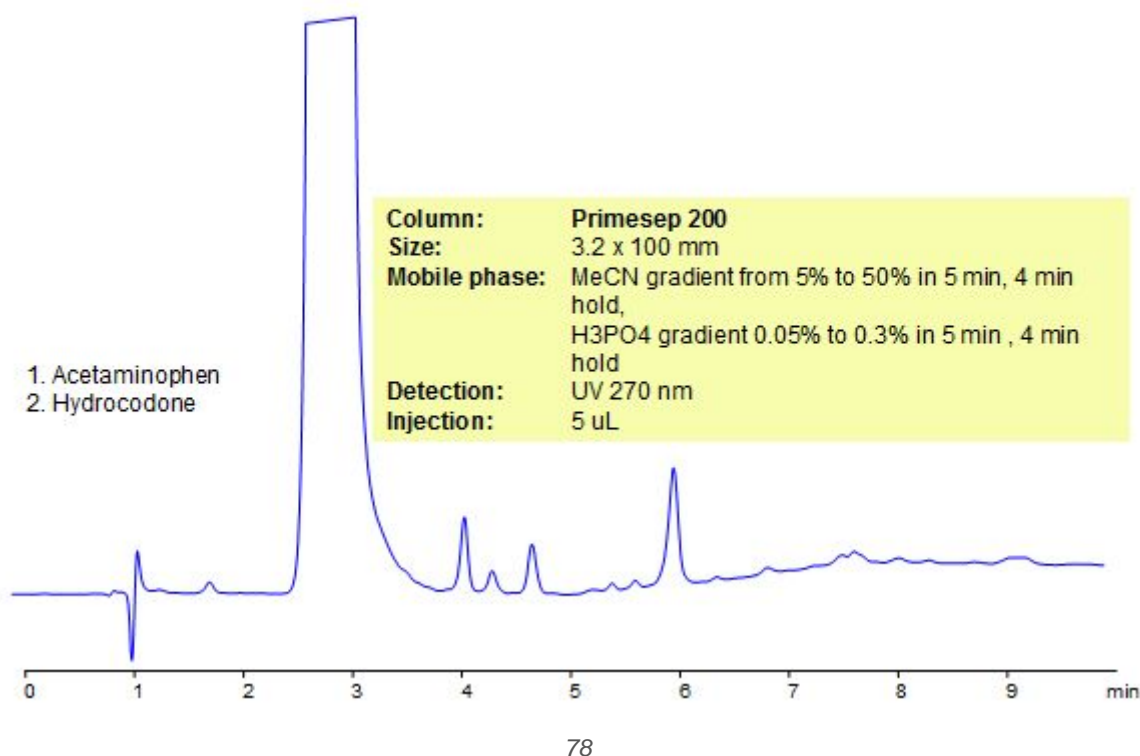

Method Parameters

Mobile Phase	Gradient MeCN – 5-50%
Buffer	Gradient H3PO4 – 0.05-0.3%, 5 min, 4 min hold
Flow Rate	0.5 ml/min
Detection	UV, 270 nm
Class of Compounds	Drug, Analgetic, Acid, Hydrophilic, Ionizable,
Analyzing Compounds	Acetaminophen, Hydrocodone

HPLC Column Used

Primesep 200, 3.2x100 mm, 5 µm, 100A

[Order this column at hplc-shop.de](http://hplc-shop.de) →