


HPLC Determination of Prochloraz on Newcrom B Column

<https://sielc.com/hplc-determination-of-prochloraz-on-newcrom-b-column>

Chromatogram

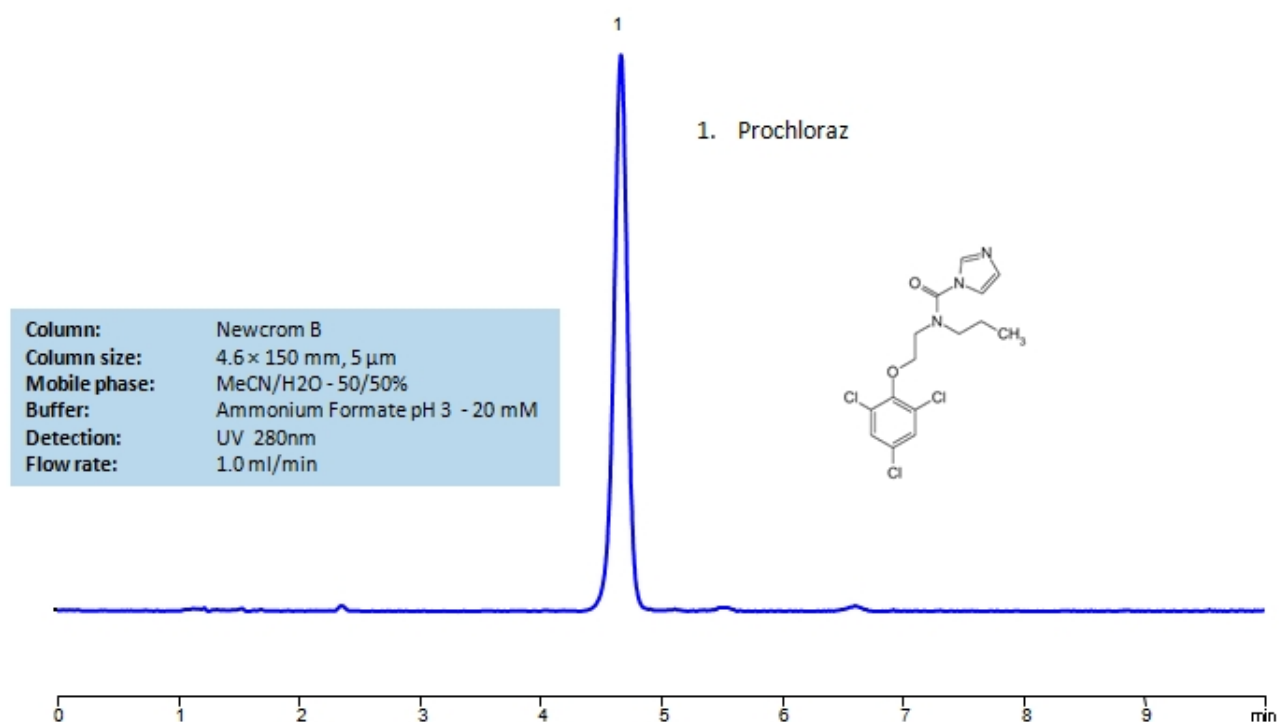

HPLC Determination of Prochloraz on Newcrom B Column_1215

Description

· High Performance Liquid Chromatography (HPLC) Method for Analysis of Prochloraz .

Prochloraz is an imidazole fungicide with the chemical formula C₁₅ H₁₆ Cl₃ N₃ O₂ . It is used widely across Europe, Asia, South America, and Australia, but is not registered for use in the United States. It works as an inhibitor of the enzyme lanosterol 14α-demethylase.

Prochloraz can be retained and analyzed using the Newcrom B stationary phase column. The analysis utilizes an isocratic method with a simple mobile phase consisting of water and acetonitrile (MeCN) with an ammonium formate buffer. Detection is performed using UV.

Method Parameters

Mobile Phase	MeCN/H ₂ O – 50/50%
Buffer	AmFm pH 3.0 – 20 mM
Flow Rate	1.0 ml/min
Detection	UV, 280 nm
Class of Compounds	Hydrophobic, Herbicide, Pesticide

HPLC Column Used

Newcrom B, 4.6 x 150 mm, 5 µm, 100 Å, dual ended[Order this column at hplc-shop.de →](https://www.hplc-shop.de)