


HPLC Method for Analysis of Trimethylglycine (Betaine) on Primesep 100 Column

<https://sielc.com/hplc-method-for-analysis-trimethylglycine>

Chromatogram

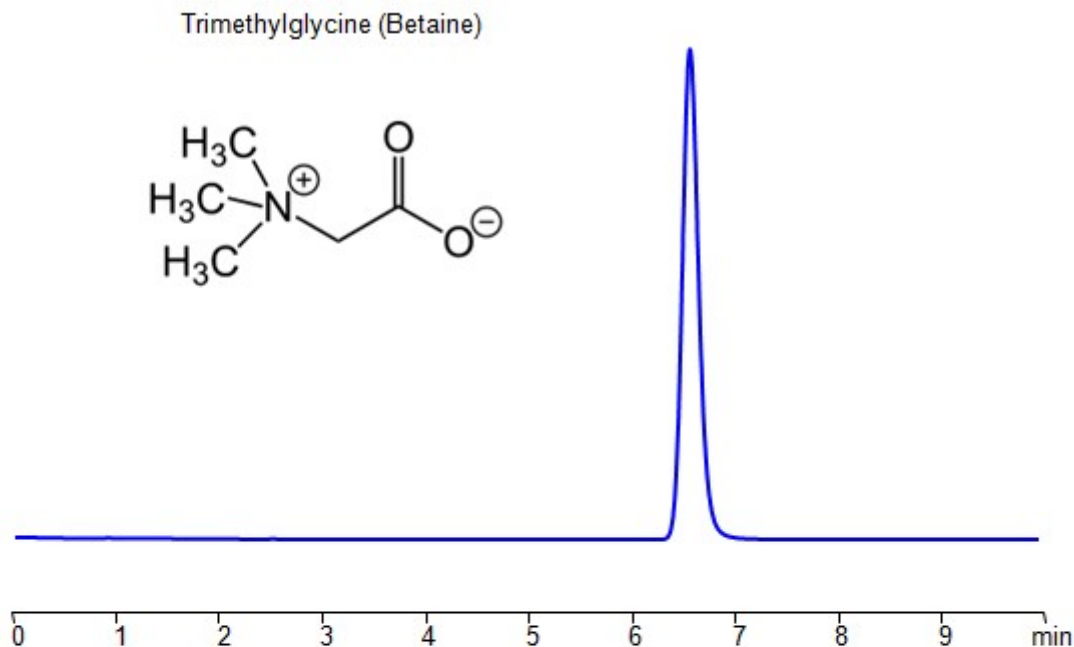

Column:	Primesep 100
Column size:	4.6 × 150 mm, 5 µm
Column part number:	100-46.150.0510
Mobile phase:	MeCN/H ₂ O – 10/90%
Buffer:	H ₂ SO ₄ - 0.1%
Flow rate:	1.0 mL/min
Detection:	UV 200 nm
Injection volume:	5.0 µL
Sample:	1.0 mg/ml
LOD:	0.8 ppm

Description

· HPLC Method for Analysis of Trimethylglycine (Betaine) on Primesep 100 Column by SIELC Technologies

Trimethylglycine, also known as betaine, is a compound derived from the amino acid glycine. It is found in various foods, including beets, spinach, whole grains, and seafood. Trimethylglycine has been studied for its potential health benefits, and it is commonly used as a dietary supplement.

Methylation Support: Trimethylglycine plays a role in a process called methylation, which is important for various biochemical reactions in the body. Methylation is involved in the synthesis of neurotransmitters, DNA, and proteins. **Heart Health:** Some studies suggest that trimethylglycine may have potential benefits for heart health. It is believed to help lower levels of

homocysteine, an amino acid linked to an increased risk of heart disease when present in high concentrations.

Exercise Performance: There is some evidence suggesting that trimethylglycine supplementation may enhance exercise performance. It is thought to improve oxygen utilization and reduce fatigue during physical activity.

Liver Health: Trimethylglycine has been investigated for its potential to support liver health. It may help protect the liver from certain forms of damage.

Anti-Inflammatory Properties: Some studies propose that trimethylglycine may have anti-inflammatory effects, which could be beneficial in various health conditions.

Trimethylglycine be retained and analyzed using a Primesep 100 mixed-mode stationary phase column. The analysis employs an isocratic method with a simple mobile phase comprising water, acetonitrile (MeCN), and sulfuric acid as a buffer. This method allows for detection using UV at 200 nm

Method Parameters

Mobile Phase	MeCN/H ₂ O – 10/90%
Buffer	H ₂ SO ₄ -0.2%
Flow Rate	1.0 ml/min
Detection	UV 200 nm
Samples	1 mg/mL
Injection volume	5 µl
LOD*	0.8 ppm
Class of Compounds	Aromatic amines, Ureas
Analyzing Compounds	Trimethylglycine (Betaine)

HPLC Column Used

Primesep 100, 4.6 x 150 mm, 5 µm, 100 A, dual ended

[Order this column at hplc-shop.de](http://hplc-shop.de) →