


HPLC Separation of Adenosine Mono-, Di- and Triphosphate on Newcrom B column

<https://sielc.com/hplc-separation-of-adenosine-mono-di-and-triphosphate-on-newcrom-b>

Chromatogram

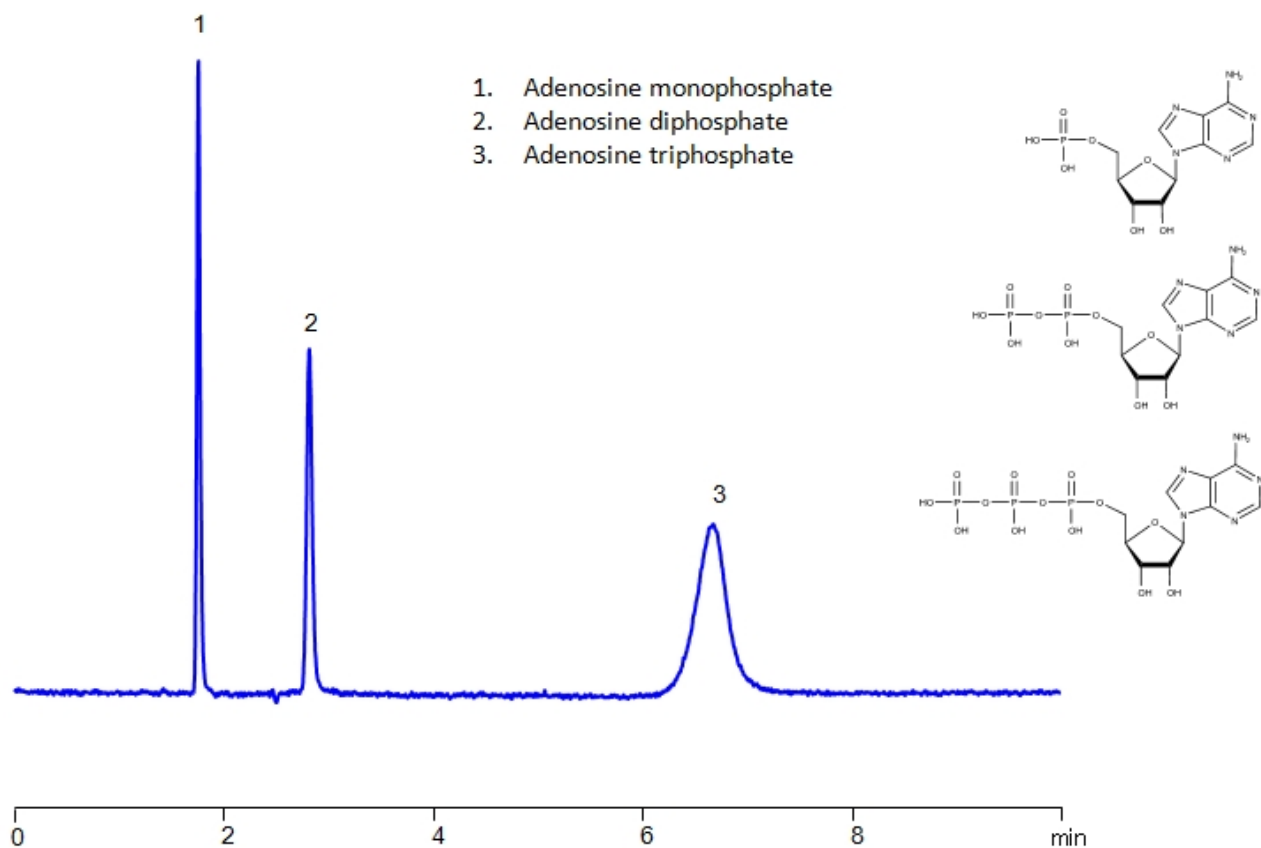

Column:	Newcrom B
Column size:	4.6 × 150 mm, 3 μm
Column part number:	NB-46.150.0310
Mobile phase:	MeCN/H ₂ O – 10/90%
Buffer:	H ₃ PO ₄ – 1.0 %
Flow rate:	1.0 mL/min
Detection:	UV 262 nm

Description

High Performance Liquid Chromatography (HPLC) Method for Analysis of Adenosine Diphosphate , Adenosine Monophosphate , Adenosine Triphosphate .

Adenosine Triphosphate (ATP) is a nucleotide with the chemical formula C₁₀ H₁₆ N₅ O₁₃ P₃. It is a primary energy carrier of the cell and is crucial to cellular energy metabolism. ATP is generated in the mitochondria and is used in numerous various cellular processes. It releases energy when a phosphate linkage breaks from it, creating adenosine diphosphate.

Adenosine Diphosphate (ADP) , also known as adenosine pyrophosphate, is a compound with the formula C₁₀ H₁₅ N₅ O₁₀ P₂ . It is a precursor in the synthesis of DNA and RNA. Medications that block ADP receptors on platelets prevent blood clots in conditions like heart attacks and strokes.

Adenosine Monophosphate (AMP) is a nucleotide with the chemical formula $C_{10}H_{14}N_5O_7P$. Due to being a byproduct of ATP and ADP, it can be reused by the body for energy at higher forms. AMP can be converted into cyclic adenosine monophosphate (cAMP) to become a second messenger in the body, relaying signals between cells.

Adenosine Diphosphate, Adenosine Monophosphate, Adenosine Triphosphate can be retained and analyzed using the Newcrom B stationary phase column. The analysis utilizes an isocratic method with a simple mobile phase consisting of water and acetonitrile (MeCN) with a Sulfuric Acid buffer. Detection is performed using UV.

Method Parameters

Mobile Phase	MeCN/H ₂ O – 10/90%
Buffer	H ₂ SO ₄ – 1.0%
Flow Rate	1.0 ml/min
Detection	UV 262nm
Class of Compounds	Acid, Hydrophilic
Analyzing Compounds	Adenosine Diphosphate, Adenosine Monophosphate, Adenosine Triphosphate

HPLC Column Used

Newcrom B, 4.6 x 150 mm, 3 µm, 100 Å, dual ended

[Order this column at hplc-shop.de →](http://hplc-shop.de)