


HPLC Method for Analysis of Tetrahydrozoline Hydrochloride in Redness Relief Eye Drops on Newcrom R1 column

https://sielc.com/hplc_method-for-analysis-tetrahydrozoline

Chromatogram

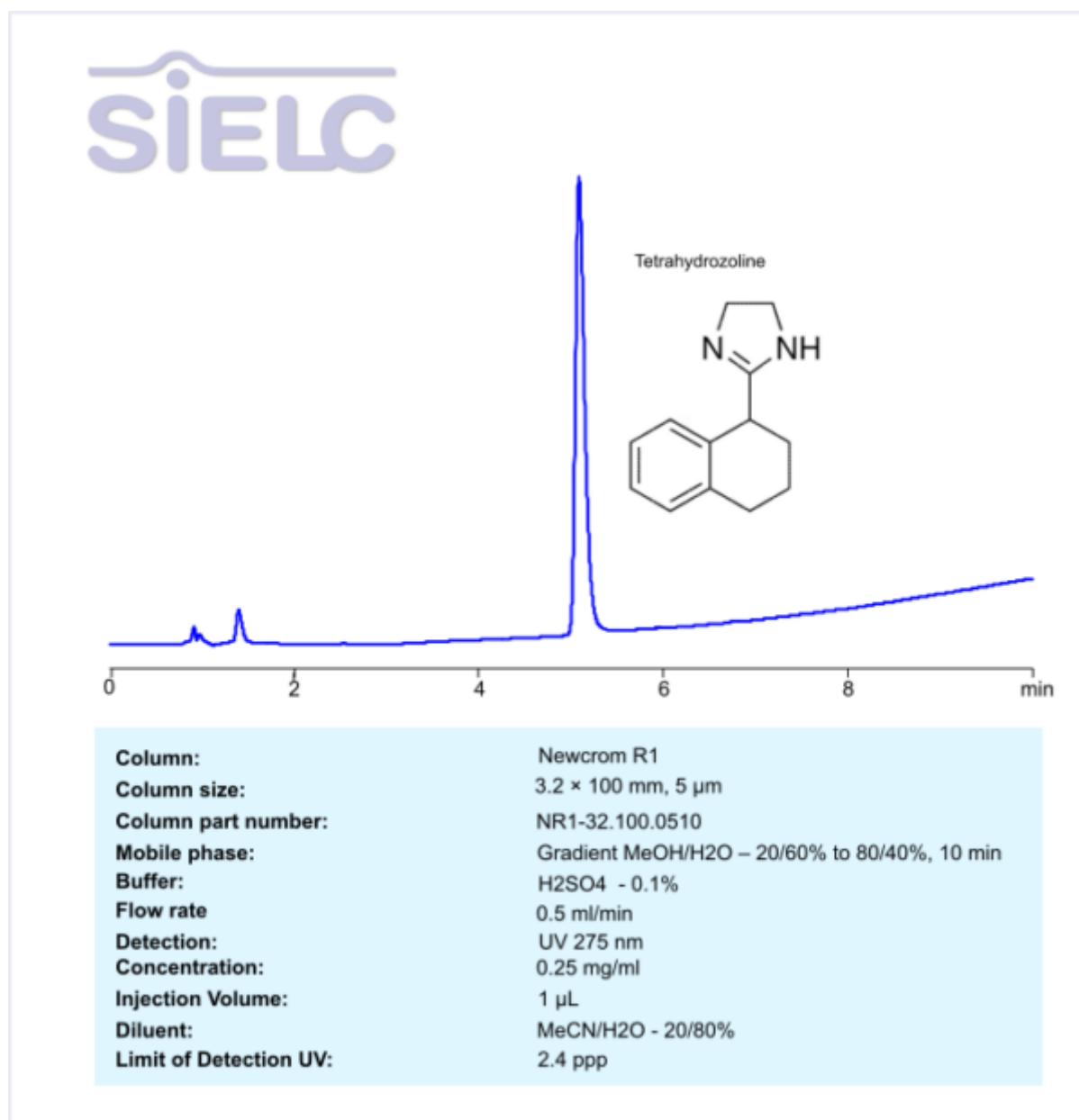

Description

· High Performance Liquid Chromatography (HPLC) Method for Analysis of Tetrahydrozoline

Tetrahydrozoline is a medication primarily used as a topical decongestant . It is commonly found in over-the-counter eye drops (e.g., Visine) and nasal sprays to relieve redness in the eyes or nasal congestion.

Tetrahydrozoline can be retained and analyzed using the Newcrom R1 stationary phase column. The analysis utilizes a gradient method with a simple mobile phase consisting of water, methanol (MeOH), and sulfuric acid. Detection is performed using UV at 275 nm.

Method Parameters

Mobile Phase	Gradient MeOH – 20-60%, 10 min
Buffer	H2SO4 – 0.1%
Flow Rate	0.5 ml/min
Detection	UV 275 nm
LOD	2.4 ppb
Class of Compounds	Drug
Analyzing Compounds	Tetrahydrozoline

HPLC Column Used

Newcrom R1, 3.2 x 100 mm, 5 µm, 100 A, dual ended

[Order this column at hplc-shop.de →](https://www.hplc-shop.de)