


# UV-Vis Spectrum of Glycine

<https://sielc.com/uv-vis-spectrum-of-glycine>

## Chromatogram

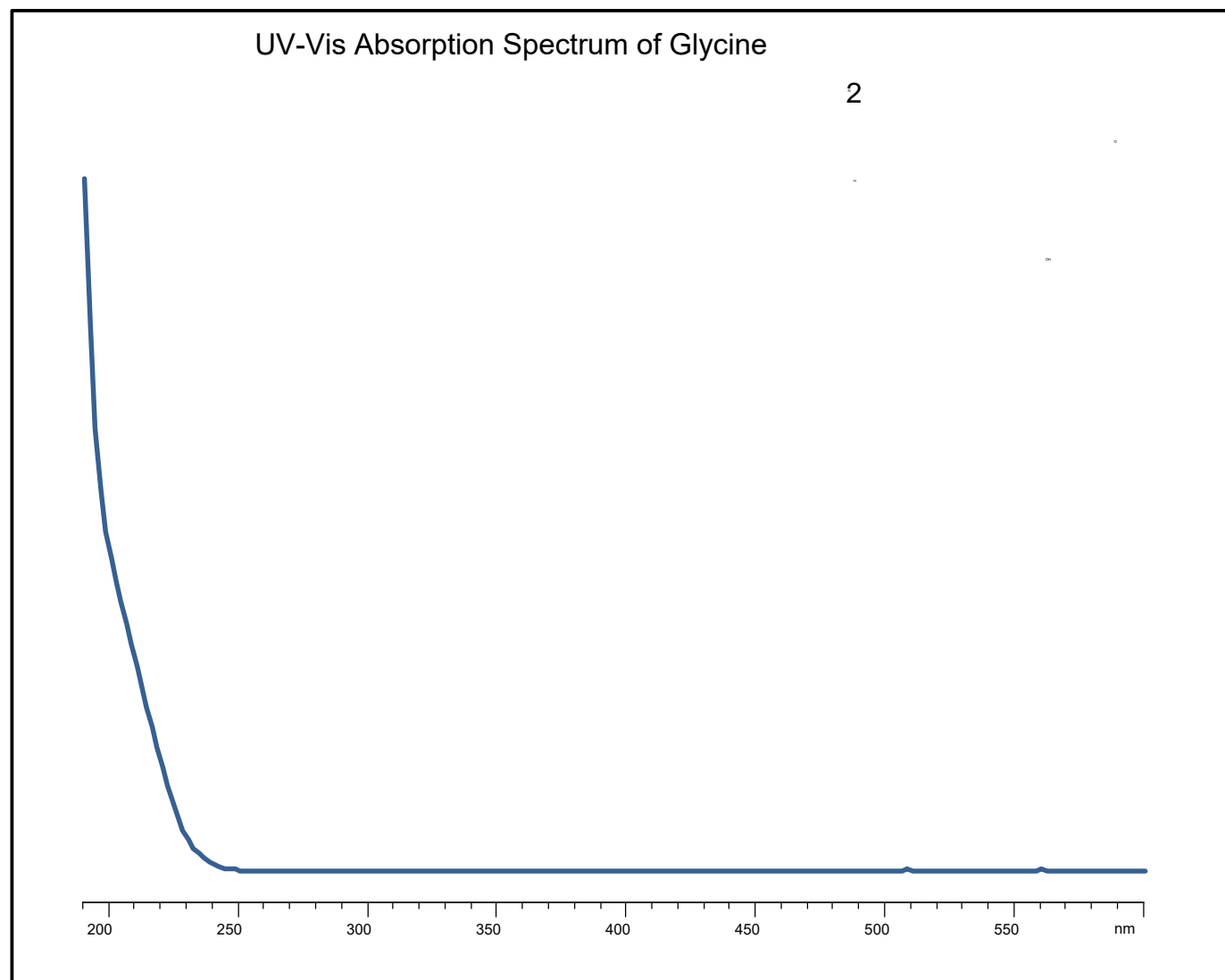

*UV-Vis Spectrum of Glycine. Absorption starts high at 200 nm but quickly falls off and bottoms out around 240 nm.*

## Description

· Access the UV-Vis Spectrum SIELC Library

For optimal results in HPLC analysis, it is recommended to measure absorbance at a wavelength that matches the absorption maximum of the compound(s) being analyzed. The UV spectrum shown can assist in selecting an appropriate wavelength for your analysis. Please note that certain mobile phases and buffers may block wavelengths below 230 nm, rendering absorbance measurement at these wavelengths ineffective. If detection below 230 nm is required, it is recommended to use acetonitrile and water as low UV-transparent mobile phases, with phosphoric acid and its salts, sulfuric acid, and TFA as buffers.